



MEDICAL TERMINAL UNITS **GCE MEDIUNIT-BSI TYPE**

EDITION 1/2014

MEDICAL TERMINAL UNITS

Medical terminal units provide quick and easy connection of hospital ward gas equipment to the hospital gas source.

GCE complies with ISO 7396 and national installation standards with secure products where every product is fully tested in production. Our Medical gas outlets are in accordance with ISO EN 9170-1 international standards.

KEY BENEFITS

- Wall housing is compatible with all GCE MediUnit standards like DIN, BSI, SS, CZ
- After 10 years it is possible to upgrade the units with a special upgrade pack
- Small size
- All functional components from brass
- Simple installation
- Easy to clean
- Fast connection and disconnection
- Attractive design
- Recessed, Exposed and Bed head installation versions
- Complies with colour coding and description by standard
- In accordance with ISO 9170-1
- Ultra safe product where every piece is fully tested in production




Recessed version



Exposed version

MEDICAL TERMINAL UNIT - BSI

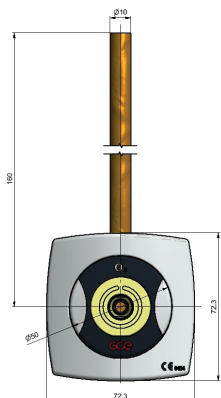

TECHNICAL DATA

Gases:	O ₂ ; O ₂ /N ₂ O; N ₂ O; Air; Air-800; VAC
Dimensions:	
Height	73 mm
Width	73 mm
Depth	63 mm
Weight:	0,36 kg
Working pressure:	4 - 5 bar
Regulatory status:	Complies with Medical directive 93/42/EEC
	Complies with EN ISO 7396-1 (Central Gas Supply system)
	Complies with EN ISO 9170-1 (Terminal units)
	Complies with BS 5682 (BSI gas specific connections)

Art. Nr.	Description	Gas	Inlet
0732046	TU BSI - recessed	O ₂	pipe Ø10
0732047	TU BSI - recessed	Air	pipe Ø10
0732048	TU BSI - recessed	VAC	pipe Ø10
0732049	TU BSI - recessed	N ₂ O	pipe Ø10
0732050	TU BSI - recessed	O ₂ /N ₂ O	pipe Ø10
0732051	TU BSI - recessed	Air-800	pipe Ø10
0732052	TU BSI - exposed	O ₂	pipe Ø10
0732053	TU BSI - exposed	Air	pipe Ø10
0732054	TU BSI - exposed	VAC	pipe Ø10
0732055	TU BSI - exposed	N ₂ O	pipe Ø10
0732056	TU BSI - exposed	O ₂ /N ₂ O	pipe Ø10
0732057	TU BSI - exposed	Air-800	pipe Ø10



Installation plug

INSTALLATION TOOLS

MP_00345	QC installation key
MP_00324	Button remover
0732040	Installation plug (10 pcs)

MOUNTING



- Gas specific connection Probe - Terminal unit
- Colour coding marking on terminal unit



- Push probe to working position like on the picture
- In working position Probe is connected to the gas supply
- Pressure drop is in accordance with EN 9170-1



- Depress the button to remove probe

UNITED KINGDOM & IRELAND GCE Ltd

Yew Tree Way, Stone Cross Park
Golborne, Warrington
WA3 3JD
Phone: +44 (0)1942 29 29 50
Fax: +44 (0)1942 29 29 77
sales.gb@gcegroup.com
www.gcegroup.com



Gas Control Equipment